



## Schedule C Planned Maintenance

PM Schedule Sys\_C

Site: <u>Green light img</u>	Month: <u>Oct</u>
System ID#: _____	Dispatch/Job #: _____
Service Engineer: <u>AG</u>	


PM Check	Test Purpose	Pass	Fail	Comments
<b>Image Quality</b>				
<u>Alignment Light Check</u>	IQ	✓		
<u>Perform SaveInfo</u>	NA	✓		
<u>Patch Update (Service contract customers only)</u>				
<u>LVShim</u>	IQ	✓		
<u>Eddy Current Check</u>	IQ	✓		
<u>Coherent Noise</u>	IQ	✓		
<u>Signal to Noise Check</u>	IQ	✓		
<u>Spike Noise Check</u>	IQ	✓		
<u>MCR Tool Check (for 14.0M5 and earlier software only)</u> OR <u>16-Channel Switch Diagnostics (for 14.0M5A and forward software only)</u>	IQ			
<u>PM Check (Service contract customers only)</u>	IQ			
<b>Magnet Room</b>				
<u>PAC Leakage Current Test</u>	Safety			
<u>Inspect Front Cover Cable Take-Up (Saltbox style magnets only)</u>	Check			
<u>Cardiac Gating Cable</u>	Safety			
<u>Patient Blower &amp; Filter</u>	Check	✓		<u>cleaned</u>
<b>RF/System Cabinet</b>				
<u>Check/Clean Filters in RFS Cabinet</u>	Clean	✓		<u>cleaned</u>
<u>Fans &amp; Filter for SRFD</u>	Clean	✓		
<u>Fan &amp; Belt for MKS</u>	Clean	✓		
<u>Check RF Output Power</u>	Check	✓		
<b>Patient Handling</b>				
<u>Patient Table Checks</u>	Safety	✓		<u>cleaned, replaced 1 wheel</u>



Lite Patient Transport: 	Safety			
<b>Gradients</b>				
<i>Check Fluid Levels and Valve of Heat Exchanger</i>	Check			
<i>Fans &amp; Filter</i>	Clean	<input checked="" type="checkbox"/>		
Lytron System Maintenance: 				
<b>Twin Accessory Cabinet</b>				
<i>Filter Replacement</i>	Replace			
<i>Clean Inlet Screen</i>	Clean			
<i>Solenoid Valve Replacement</i>	Replace			
<b>Computer</b>				
<i>Storelog</i>	Check	<input checked="" type="checkbox"/>		
<b>Magnet</b>				
<i>Test GE Magnet Rundown Unit (MRU)</i>	Safety			
<b>Coils</b>				
<i>Coil Shock Sensor Check</i>	Check			

Schedule C Signa Planned Maintenance.

Date Completed: 10/7/22

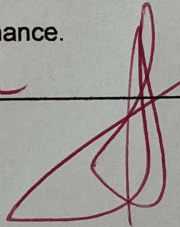

8ch torso - needs repair  
 knee coil - needs repair  
 PBS hose connector - damaged  
 Replaced 2 wheel  
 Cradle latch, cracked.  
 CBL coil - needs repair.

Back to top