



Schedule D Planned Maintenance

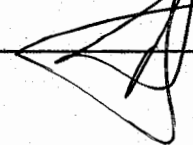
PM Schedule Sys_D

Site: <u>Green light IMG</u>	Month: <u>1/20/2013</u>
System ID#: _____	Dispatch/Job #: _____
Service Engineer: <u>A.G.</u>	


PM Check	Test Purpose	Pass	Fail	Comments
Image Quality				
<i>Alignment Light Check</i>	IQ	✓		
<i>Perform SaveInfo</i>	N/A	✓		
<i>Patch Update (Service contract customers only)</i>				
<i>LVShim</i>	IQ	✓		
<i>Eddy Current Check</i>	IQ	✓		
<i>Coherent Noise</i>	IQ	✓		
<i>Signal to Noise Check</i>	IQ	✓		
<i>Spike Noise Check</i>	IQ	✓		
<i>MCR Tool Check (for 14.0M5 and earlier software only)</i>				
OR				
<i>16-Channel Switch Diagnostics (for 14.0M5A and forward software only)</i>	IQ			
<i>PM Check (Service contract customers only)</i>	IQ			
Magnet Room				
<i>Pneumatic Patient Alert System</i>	Safety			
<i>Cardiac Gating Cable</i>	Safety			
<i>Patient Blower & Filter</i>	Check	✓		<i>Cleared</i>
RF/System Cabinet				
<i>Fans & Filter for SRFD</i>	Clean	✓		<i>cleaned</i>
<i>Fan & Filter for MKS</i>	Clean	✓		<i>cleaned</i>
<i>Check RF Output Power</i>	Check			
Patient Handling				
<i>Patient Table Checks</i>	Safety	✓		<i>Cleaned and checked</i>
<i>Lite Patient Transport:</i>	Safety			
Gradients				

CFT-150 System Maintenance:				
				
Mobile Hx System Maintenance:				
				
Contact van owner for maintenance on cooling systems in addition to Mobile Hx (LTL-4, M1085KF).				
Fans & Filter		Clean	<input checked="" type="checkbox"/>	Cleaned
Twin Accessory Cabinet				
Filter Replacement		Replace	<input type="checkbox"/>	
Computer				
Storelog		Check	<input checked="" type="checkbox"/>	
Clean Dust, Set Time		Clean	<input checked="" type="checkbox"/>	Cleaned, set time
Magnet				
Test GE Magnet Rundown Unit (MRU)		Safety	<input type="checkbox"/>	
Inspect GE Magnet Rundown Unit (MRU)		Safety	<input type="checkbox"/>	
Inspecting Cryocooler Systems		Check	<input type="checkbox"/>	
Main Disconnect Panel (MDP)				
Check battery voltage. Refer to MDP Manual. (Need certified electrician with Arc Flash PPE to perform this task.)		Check	<input type="checkbox"/>	
Coils				
Coil Shock Sensor Check		Check	<input checked="" type="checkbox"/>	

Schedule D Signa Planned Maintenance.

Date Completed: 1/20/23 

Fixed door
Needs Head coil for Q4.

 Back to top